


Genotyping Report

Strain ID	T028553	Strain Type	KO(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Yin Chen	Gene Name	Adcy6		

1. Strategy of Genotyping

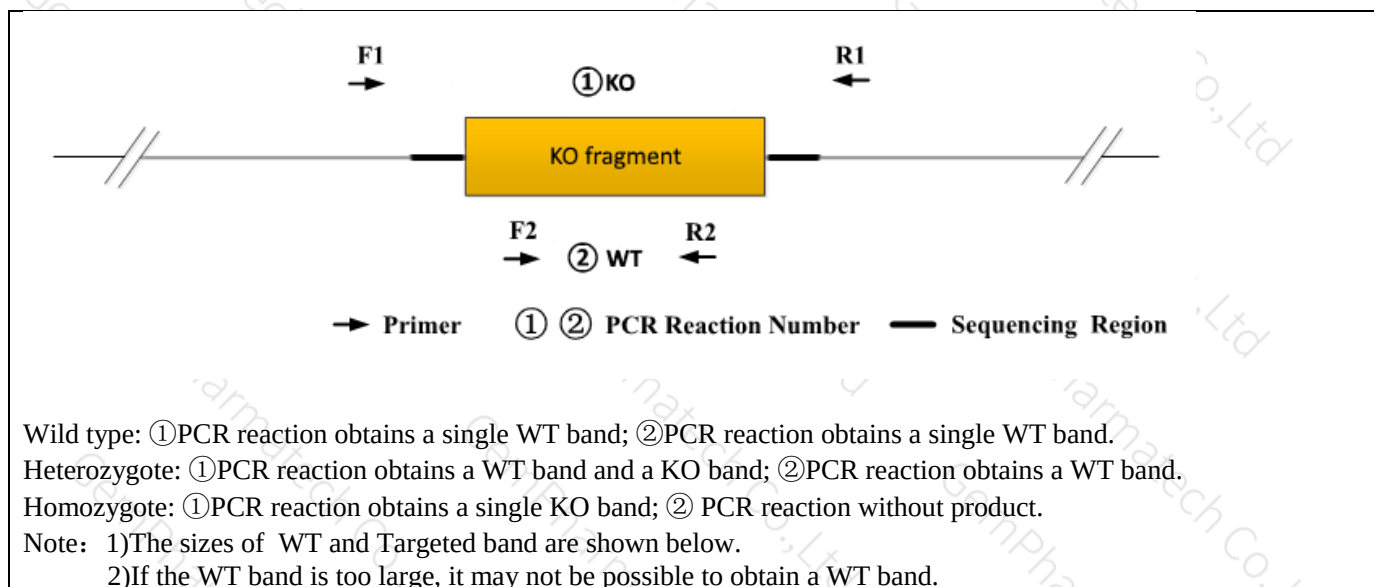


2. Primer Information

PCR No.	Primer No.	Primer Name	Sequence	Band Size
PCR①	F1	JS09486-Adcy6-5wt-tF1	GCTTGGGAATCTAAGCCAAGG	WT:8702bp KO: 302bp
	R1	JS09486-Adcy6-3wt-tR1	ACTGAAGGCTTCTCCAGGCATT	
PCR②	F2	JS19486-Adcy6-wt-tF1	AAGATCAAGACCATCGGTAGCAC	WT:369bp KO:0bp
	R2	JS19486-Adcy6-wt-tR1	CTTAGAATGTACTGGACACGTCTAG G	

3. Gel Image

ggaggctttagagatatctcctcctgtt-----8400bp-----tgggtggaaggcagtgattaataacatacg
--

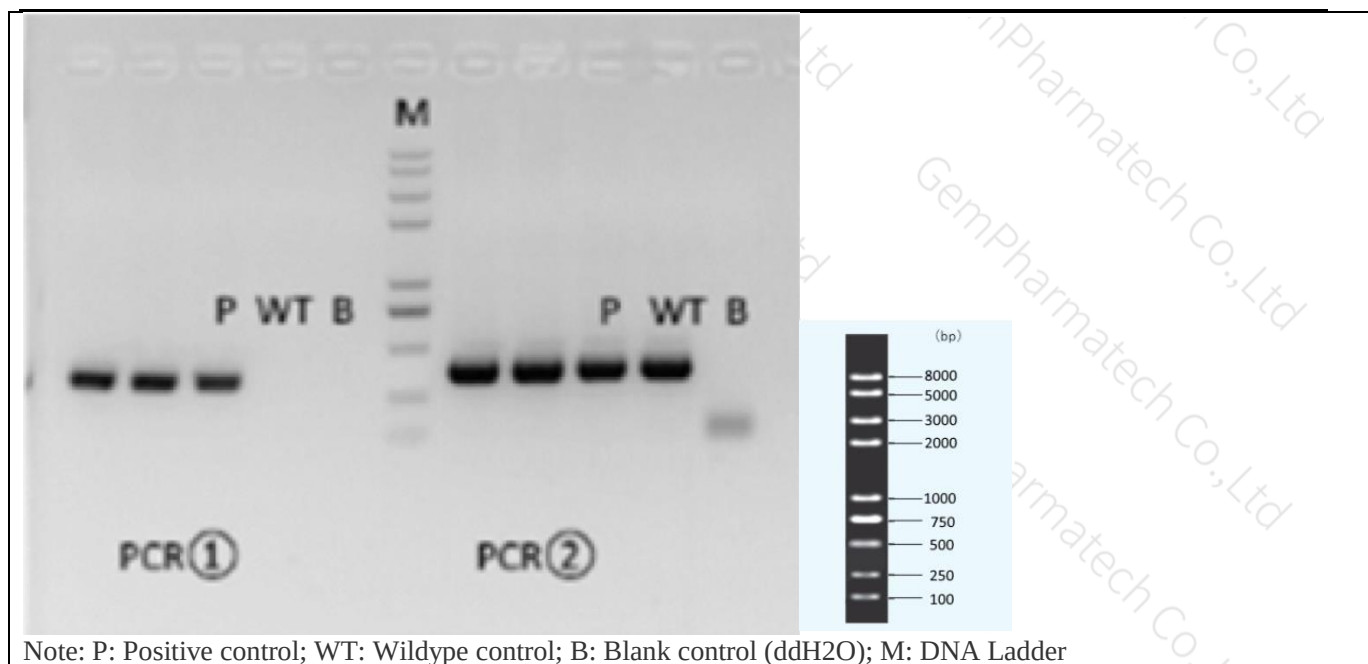

4. PCR Condition

PCR Reaction Component			
Seg.	reaction component		Volume (μl)
1	2 × Taq Master Mix , Dye Plus, (Vazyme P112-03)		12.5
2	ddH2O		9.5
3	Primer A(10pmol/μl)		1
4	Primer B(10pmol/μl)		1
5	Template(≈100ng/μl)		1
PCR program			
Seg.	Temp.	Time	Cycle
1	95℃	5min	20×
2	98℃	30s	
3	65℃ (-0.5℃/cycle)	30s	
4	72℃	45s	
5	98℃	30s	20×
6	55℃	30s	
7	72℃	45s	
8	72℃	5min	
9	10℃	hold	

Reviewer:Ting Sun

Date:2020-7-1



江苏集萃药康生物科技股份有限公司

GemPharmatech Co., Ltd