


Genotyping Report

Strain ID	T017460	Strain Type	KO(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Zifan Lin	Gene Name	Col4a3bp		

1. Strategy of Genotyping


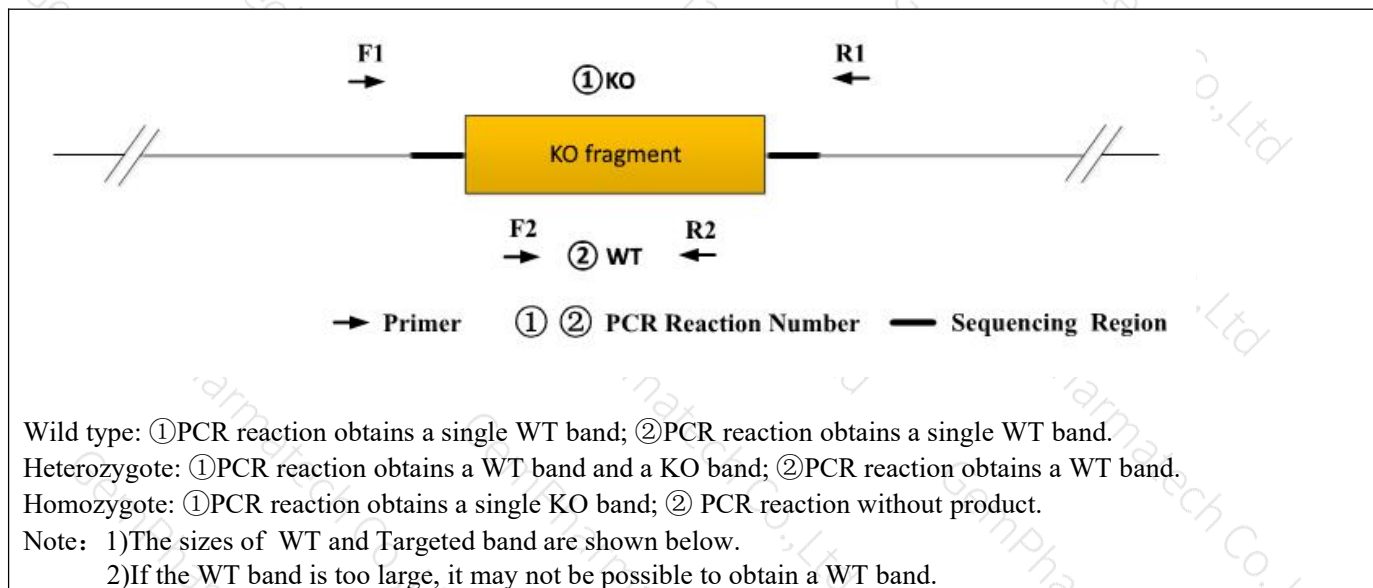


2. Primer Information

PCR No.	Primer No.	Sequence	Band Size
PCR①	T017460-F1	GTCAGGTATTGAGTGACCTGTATGACTT	WT:53037bp KO: 445bp
	T017460-R1	AGCTCTATGACTTCTCCAGTTGCC	
PCR②	T017460-F2	GGTTCCCAGATCAATGTTATCTGTCC	WT:399bp KO:0bp
	T017460-R2	GGATAACATCCTCCTGACTGATCG	

3. Gel Image

aacaggaacatagggaggaattacctcacgggag---52592bp---taagaactgggtaccactccctggaatcttttca

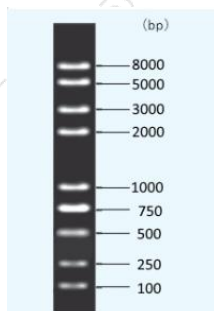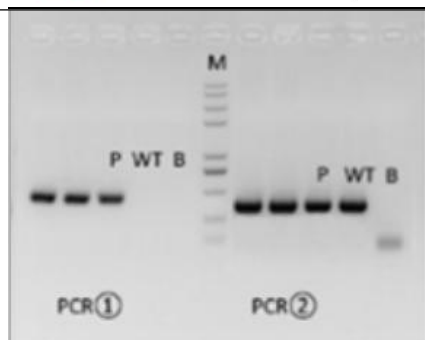

Note: P: Positive control; WT: Wildtype control; B: Blank control (ddH₂O); M: DNA Ladder

4. PCR Condition

PCR Reaction Component			
Seg.	reaction component		Volume (μl)
1	2 × Taq Master Mix , Dye Plus, (Vazyme P112-03)		12.5
2	ddH2O		9.5
3	Primer A(10pmol/μl)		1
4	Primer B(10pmol/μl)		1
5	Template(≈100ng/μl)		1
PCR program			
Seg.	Temp.	Time	Cycle
1	95℃	5min	20×
2	98℃	30s	
3	65℃ （-0.5℃/cycle）	30s	
4	72℃	45s	
5	98℃	30s	20×
6	55℃	30s	
7	72℃	45s	
8	72℃	5min	
9	10℃	hold	