


Genotyping Report

Strain ID	T003768	Strain Type	KI(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Dongdong Zhang	Gene Name	<i>Lgr5-CreERT2</i>		

1. Strategy of Genotyping

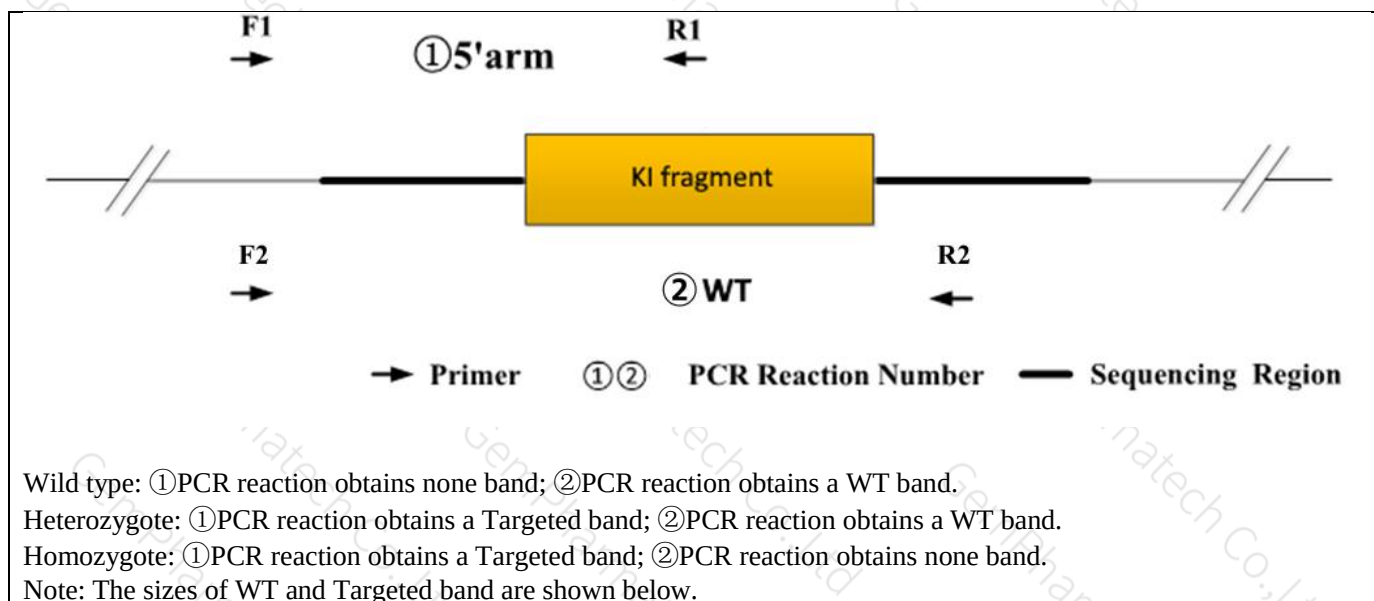


2. Primer Information

PCR No.	Primer No.	Sequence	Band Size
①5'arm	T003768-F1	CGAGTCTGCTGTCCATTAAGTC	WT:0bp Targeted:630bp
	T003768-R1	ACGCTTGAGGAGAGCCATTG	
②WT	T003768-F2	CGAGTCTGCTGTCCATTAAGTC	WT:409bp Targeted:3013bp
	T003768-R2	GTTACCTAAGTGACTTGATGGC	

3. Gel Image & Conclusion

--

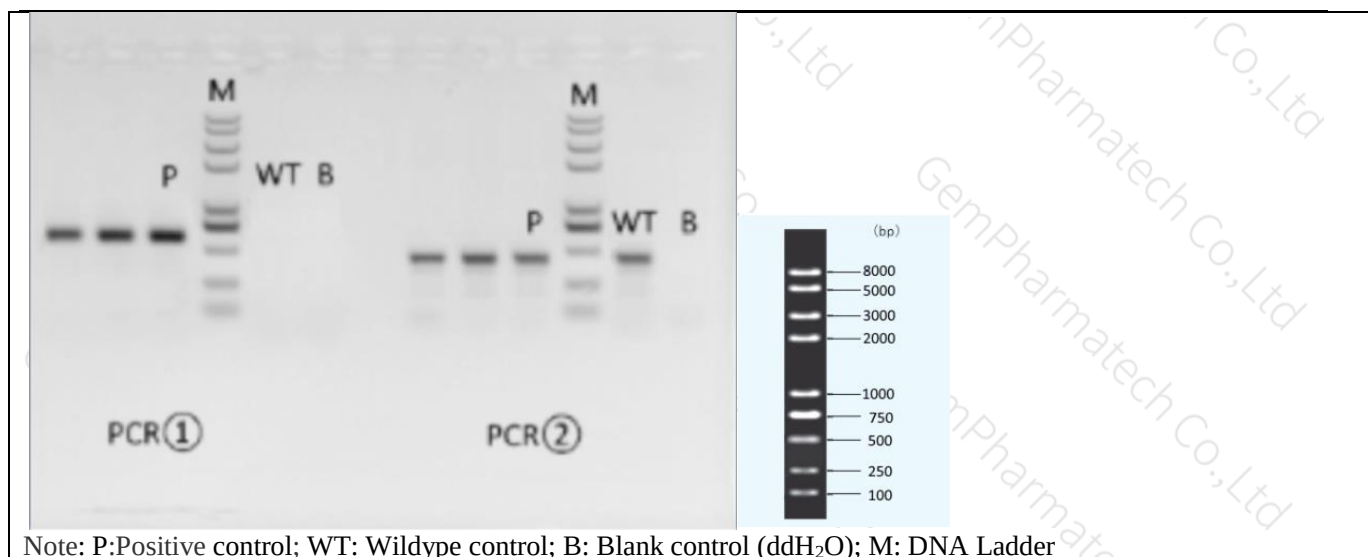

4. PCR Condition

PCR Reaction Component			
Seg.	reaction component		Volume (μl)
1	2 × Taq Master Mix , Dye Plus, (Vazyme P112-03)		12.5
2	ddH2O		9.5
3	Primer A(10pmol/μl)		1
4	Primer B(10pmol/μl)		1
5	Template(≈100ng/μl)		1
PCR program			
Seg.	Temp.	Time	Cycle
1	95℃	5min	20×
2	98℃	30s	
3	65℃ (-0.5℃/cycle)	30s	
4	72℃	45s	
5	98℃	30s	20×
6	55℃	30s	
7	72℃	45s	
8	72℃	5min	
9	10℃	hold	