

Genotyping Report

Strain ID	T009917	Strain Type	CKO(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Ya'nan Xu	Gene Name	Acss3		

1. Strategy of Genotyping



Wild type: ①PCR reaction obtains a single WT band; ②PCR reaction obtains a single WT band.

Heterozygote: ①PCR reaction obtains a WT band and a Targeted band; ②PCR reaction obtains a WT band and a Targeted band.


Homozygote: ①PCR reaction obtains a single Targeted band; ②PCR reaction obtains a single Targeted band.

Note: The sizes of WT and Targeted band are shown below.

2. Primer Information

PCR No.	Primer No.	Sequence	Band Size
①(5'arm)	T009917-F1	CATTCCTGGGAAGAGTTTGAGC	WT:263bp Targeted:368bp
	T009917-R1	GGCAAACCCCTCCCTGATAGATT	
②(3'arm)	T009917-F2	GGAAATTATTCCAGCCTTGTGG	WT:304bp Targeted:405bp
	T009917-R2	TAAGTGCGGCCAATTAGACTGC	

3. Gel Image & Conclusion

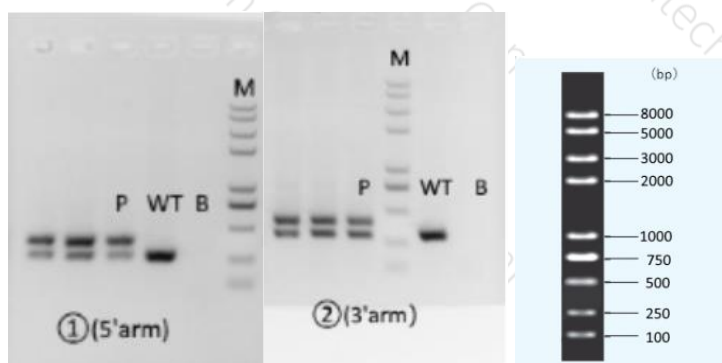

Note: P: Positive control; WT: Wildtype control; B: Blank control (ddH₂O); M: DNA Ladder

4. PCR Condition

PCR Reaction Component			
Seg.	reaction component		Volume (μl)
1	2 × Taq Master Mix , Dye Plus, (Vazyme P112-03)		12.5
2	ddH2O		9.5
3	Primer A(10pmol/μl)		1
4	Primer B(10pmol/μl)		1
5	Template(≈100ng/μl)		1
PCR program			
Seg.	Temp.	Time	Cycle
1	95℃	5min	20×
2	98℃	30s	
3	65℃（-0.5℃/cycle）	30s	
4	72℃	45s	
5	98℃	30s	20×
6	55℃	30s	
7	72℃	45s	
8	72℃	5min	
9	10℃	hold	

Reviewer:Ting Sun

Date:2021-7-23