


Genotyping Report

Strain ID	T004819	Strain Type	KI(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Dongdong Zhang	Gene Name	<i>Itgax-iCre</i>		

1. Strategy of Genotyping

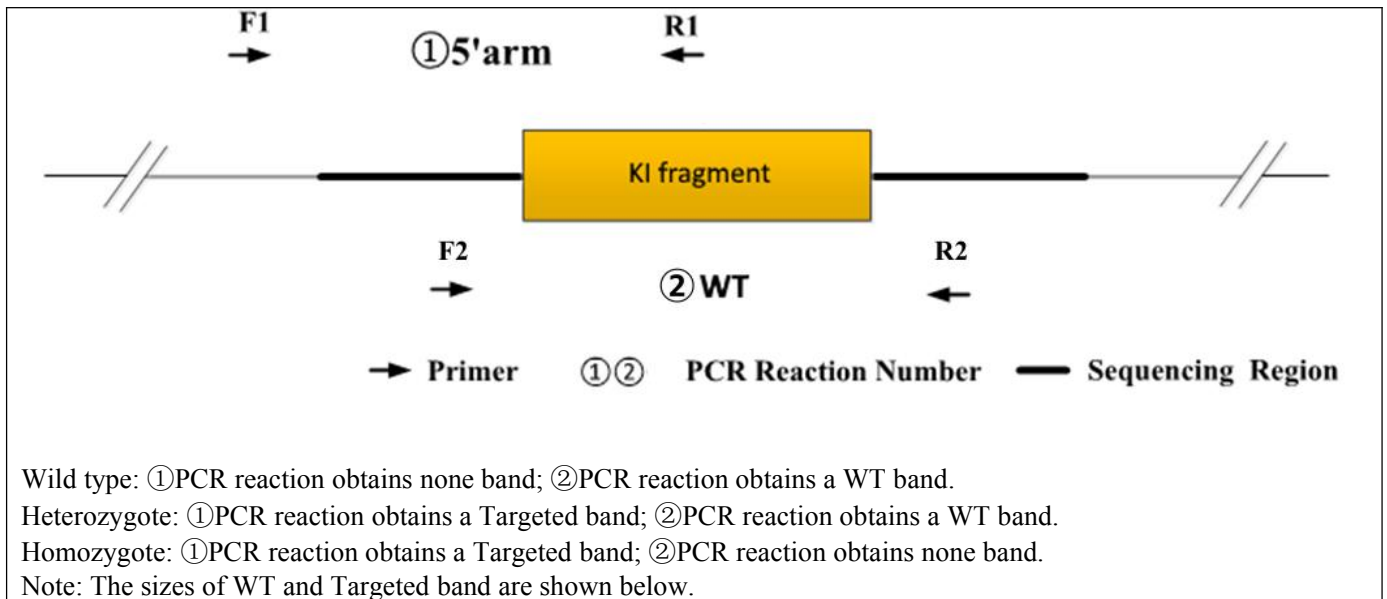


2. Primer Information

PCR No.	Primer No.	Primer Name	Sequence	Band Size
① 5' arm	F1	T004819-F1	GGGCAGTCTGGTACTTCCAAGCT	WT:0bp
	R1	T004819-R1	GCTAAGCACTGCACAGAGGATAGC	Targeted:368bp
② WT	F2	T004819-F2	CAGCAAAACCTGGCTGTGGATC	WT:412bp
	R2	T004819-R2	ATGAGCCACCATGTGGGTGTC	Targeted:0bp

3. Gel Image & Conclusion

--

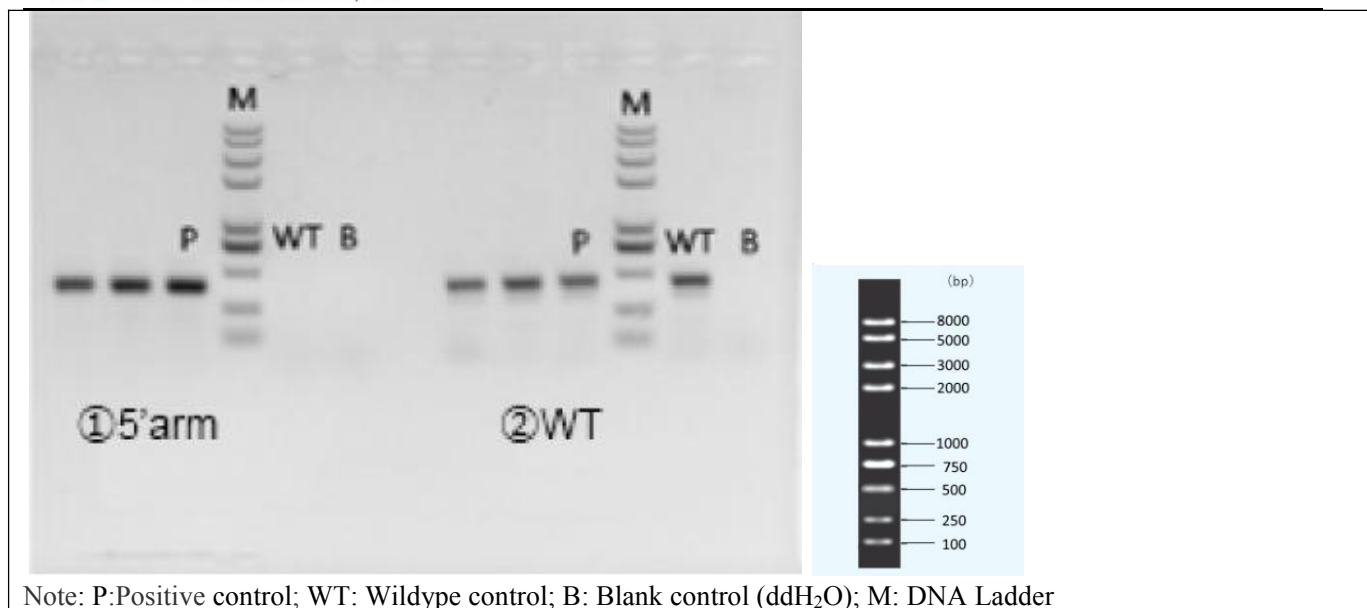

4. PCR Condition

PCR Reaction Component			
Seg.	reaction component		Volume (μl)
1	2 × Taq Master Mix , Dye Plus, (Vazyme P112-03)		12.5
2	ddH2O		9.5
3	Primer A(10pmol/μl)		1
4	Primer B(10pmol/μl)		1
5	Template(≈100ng/μl)		1
PCR program			
Seg.	Temp.	Time	Cycle
1	95℃	5min	20×
2	98℃	30s	
3	65℃ （-0.5℃/cycle）	30s	
4	72℃	45s	
5	98℃	30s	20×
6	55℃	30s	
7	72℃	45s	
8	72℃	5min	
9	10℃	hold	