


Genotyping Report

Strain ID	T006187	Strain Type	KI(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Dongdong Zhang	Gene Name	<i>Cspg4-CreERT2</i>		

1. Strategy of Genotyping

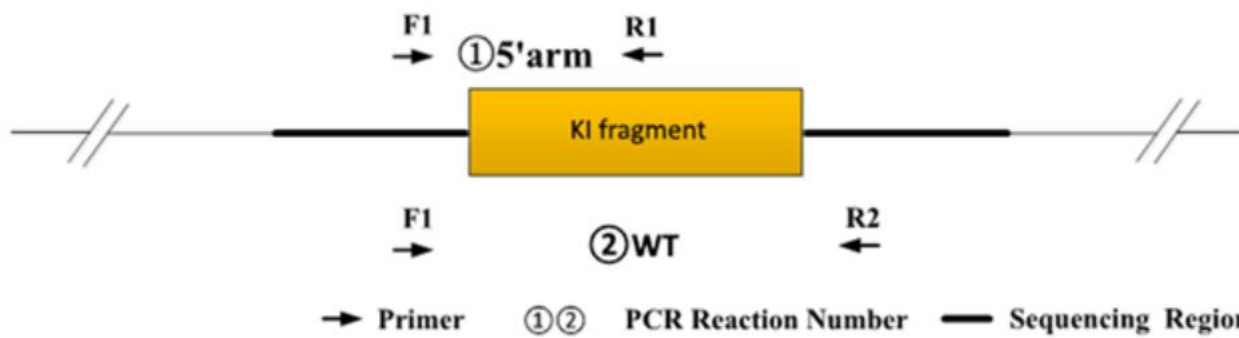


Wild type: ①PCR reaction obtains none band; ②PCR reaction obtains a WT band.
Heterozygote: ①PCR reaction obtains a Targeted band; ②PCR reaction obtains a WT band.
Homozygote: ①PCR reaction obtains a Targeted band; ②PCR reaction obtains none band.
Note: The sizes of WT and Targeted band are shown below. For ②PCR reaction, because the WT band is much smaller than the target band, it is likely to produce dominant amplification, the reaction is only used to judge whether there is a WT allele.

2. Primer Information

PCR No.	Primer No.	Primer Name	Sequence	Band Size
① 5'arm	F1	T006187-F1A	AAATCTAAGCGCGGGTCTGGC	WT:0bp Targeted:506bp
	R1	T006187-R1A	CGGACCGACGATGAAGCATGT	
② WT	F1	T006187-F1A	AAATCTAAGCGCGGGTCTGGC	WT:352bp Targeted:2386bp
	R2	T006187-R2	GGACCATGAGTGCAGTCCCCATA	

3. Gel Image & Conclusion

--

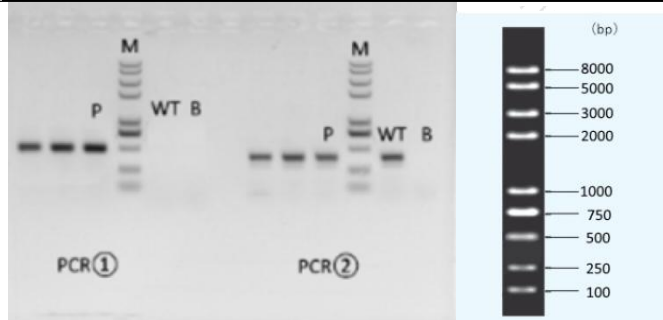

Note: P:Positive control; WT: Wildtype control; B: Blank control (ddH₂O); M: DNA Ladder

4. PCR Condition

PCR Reaction Component			
Seg.	reaction component		Volume (μl)
1	2 × Taq Master Mix , Dye Plus, (Vazyme P112-03)		12.5
2	ddH2O		9.5
3	Primer A(10pmol/μl)		1
4	Primer B(10pmol/μl)		1
5	Template(≈100ng/μl)		1
PCR program			
Seg.	Temp.	Time	Cycle
1	95℃	5min	20×
2	98℃	30s	
3	65℃ (-0.5℃/cycle)	30s	
4	72℃	45s	
5	98℃	30s	20×
6	55℃	30s	
7	72℃	45s	
8	72℃	5min	
9	10℃	hold	