


Genotyping Report

Strain ID	T005674	Strain Type	KI(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Dongdong Zhang	Gene Name	<i>Ncr1-iCre-P2A</i>		

1. Strategy of Genotyping

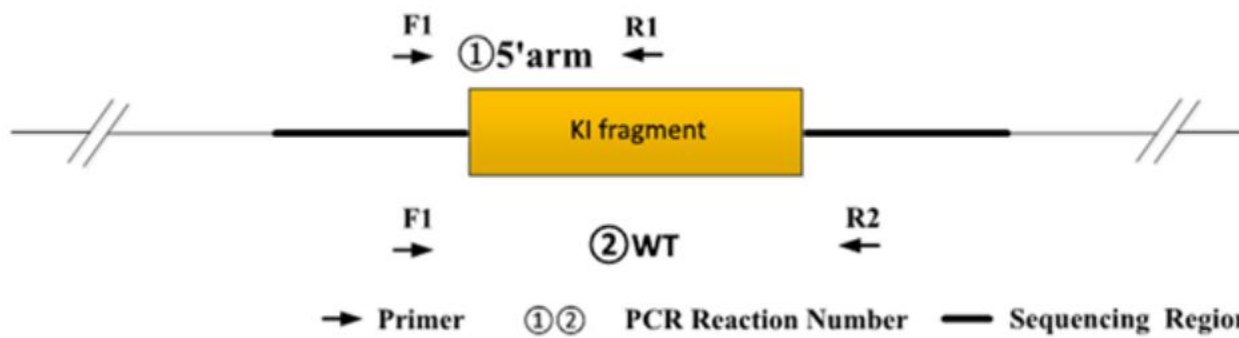


Wild type: ①PCR reaction obtains none band; ②PCR reaction obtains a WT band.
Heterozygote: ①PCR reaction obtains a Targeted band; ②PCR reaction obtains a WT band.
Homozygote: ①PCR reaction obtains a Targeted band; ②PCR reaction obtains none band.
Note: The sizes of WT and Targeted band are shown below.

2. Primer Information

PCR No.	Primer No.	Primer Name	Sequence	Band Size
①5'arm Alternative	F1	T005674-F1A	CACTGGACTCAATCCACACACATAG	WT:0bp Targeted:534bp
	R1	T005674-R1A	AGGAAGGCCAGGTTCTCTGATGT	
②WT Alternative	F1	T005674-F1A	CACTGGACTCAATCCACACACATAG	WT:444bp Targeted:1551bp
	R2	T005674-R2A	CGAGTTTCCATTGTGACCATGAT	

3. Gel Image & Conclusion

--

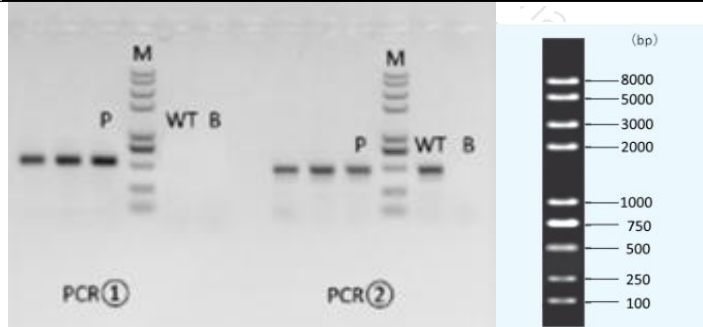

Note: P:Positive control; WT: Wildtype control; B: Blank control (ddH₂O); M: DNA Ladder

4. PCR Condition

PCR Reaction Component			
Seg.	reaction component		Volume (μl)
1	2 × Taq Master Mix , Dye Plus, (Vazyme P112-03)		12.5
2	ddH2O		9.5
3	Primer A(10pmol/μl)		1
4	Primer B(10pmol/μl)		1
5	Template(≈100ng/μl)		1
PCR program			
Seg.	Temp.	Time	Cycle
1	95℃	5min	20×
2	98℃	30s	
3	65℃ （-0.5℃/cycle）	30s	
4	72℃	45s	
5	98℃	30s	20×
6	55℃	30s	
7	72℃	45s	
8	72℃	5min	
9	10℃	hold	